

Basic Theory

Liquid Chromatography MP → Liquid
 ▶ Actively participate in equilibrium process

SP → Quasi/porous solid
 ▶ Most common do reverse-phase chromatography
 SP
 ▶ Film thickness → Very small (monolayer)

$$D_m \sim 10D_s$$

$$\text{▶ } B/U \sim 0$$

$$df \sim 0$$

$$\text{▶ } C_s U \sim 0$$

Always carried-out in packed columns

Advantage More versatile than GC

Adaptable to needs

Disadvantage Much less efficient than GC

Diffusion coefficient of analyte is orders of magnitude smaller than in GC
 ▶ Bounced into other molecules → Diffusion rate is small

In liquid phase (not gas)

Improve efficiency Use small particle → Narrow range velocity

- ▶ Dependent on particle diameter (d_p)
- ▶ Smaller particle size = smaller plate height at any given velocity

Compensate for travel distance of analyte to reach surface of SP


- ▶ Minimize the space between the particles that the analyte have to diffuse across

Analytical LC → $< 5\mu\text{m}$

HPLC → $> 3\mu\text{m}$

UHPLC → $< 3\mu\text{m}$

Theory Equations

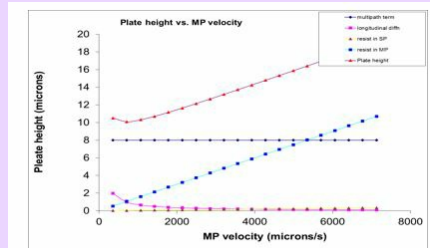


$$A = 2\lambda d_p \quad \frac{B}{U} = \frac{2\gamma D_m}{U} \approx 0$$

$$C_s U = \frac{f_s(k')d_f^2}{D_s} U \approx 0$$

$$C_m U = \frac{f_m(k')d_p^2}{D_m} U$$

LC Velocity Range

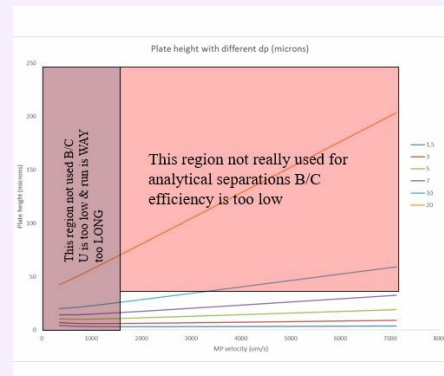

Column

Packed Column Packed full of particles

Column

Usually composed of stainless steel, etc.

Length

- ▶ 2-20cm(analytical)
- ▶ Large column → Use for preparative scale
- ▶ Small column → Packaged inside a capillary (75-100um diameter) → Couple efficiently to MS



By shaylannxd

Not published yet.

Last updated 20th June, 2022.

Page 1 of 13.

Sponsored by [ApolloPad.com](https://apollopad.com)

Everyone has a novel in them. Finish Yours!

Yours!

<https://apollopad.com>

Column (cont)

Efficiency (N)

- ▶ 3000-20000
- ▶ Dramatically lower than GC
- ▶ Rule of thumb (GUESTIMATE ONLY): $N \sim 3500 * L(\text{cm}) / d_p(\text{um})$

Sample Capacity


- ▶ Depending on size of column and packing
- ▶ ~5mg/g \rightarrow C18/silica
- ▶ ~10mg \rightarrow "vanilla" column

Thickness of SP

- ▶ ~ 1-2nm

Resistance to mass transfer in MP and multipath terms dominate

HPLC Column

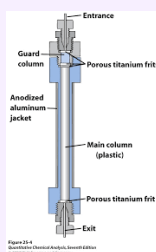

HPLC System

MP reservoirs Stores MP in inert glass bottles
 ▶ Platic coated Pyrex bottle (common) \rightarrow \$300/pc

Degas solvent

- ▶ Add element for filtering/degas
- ▶ Minimize amount of oxygen dissolved into MP \rightarrow Oxygen reactive in high pressure (increase oxidation of analyte)
- ▶ Small bubbles can form \rightarrow Result intensive undesirable peaks

Connected to a computer(pump)

- ▶ Control mixing value to produces desired MP mixture

HPLC System (cont)

Analytical Column Wide variety
 ▶ Diameter \rightarrow ~0.5cm (general)
 ▶ Length \rightarrow 10-20cm (general)

Injector Manual
 ▶ Syringe with sample
 ▶ Inject needle into port and release
 ▶ Liquid flow into loop (at atmospheric pressure)
 ▶ When rotate lever to 60 degree \rightarrow Rearrange injector (set of valves)
 ▶ Switches loop into flow path = swept down into analytical column

Operates at very high pressure

- ▶ If inject sample into septum \rightarrow Shatter syringe

Autosampler

- ▶ Prepare in vials \rightarrow Tightly seal for sample to not evaporate
- ▶ Program computer
- ▶ Runs separation overnight
- ▶ ~100 samples

Detector Record data and integrate peak area
 ▶ Quantitation

High-Pressure Pump Direct MP through system

Use high pressure

- ▶ Analytical column is filled with fine particles

Dynamic Mixer Needed to blend the different fluids (MP)

Provides the correct percentages of fluids dynamically as the separation goes



By shaylannxd

Not published yet.
 Last updated 20th June, 2022.
 Page 2 of 13.

Sponsored by **ApolloPad.com**
 Everyone has a novel in them. Finish Yours!

HPLC System (cont)

Guard Column Avoid column killers

- ▶ Species that strongly store in SP → Never eluted
- ▶ Contaminated column
- ▶ Can change separation → Destroy separation in terms of its analytical quality

Very small column

- ▶ Contains same type of SP as analytical column

Contamination is trapped inside

- ▶ Periodically replace cartridges → Preserves analytical column

Optional

Narrow Bore Tubing Tubes that connects components

Has to be rated for HPLC

- ▶ Can handle high pressure
- ▶ Has to be narrow bore → Don't want MP to mixing + dilute sample peaks

Use short length as much as possible

If fitting not installed correctly

- ▶ Can result in dead volume → Analyte that gets trapped in dead volume = gets broaden

Thermostat Oven Constant temperature


- ▶ 30-45 degrees
- ▶ For a given MP → Equilibrium is constant

Fraction Collector Robots that periodically move tubes of the eluted species from detector

Collect in vials

- ▶ Sophisticated → Deposit 1 peak per vial
- ▶ Less sophisticated → Periodically move from one vial to the next

HPLC System Diagram

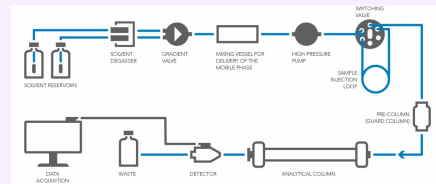

Stationary Phase

Control Retention Control retention → Control distribution constant(K)

Control by:

- ▶ Adjust type of MP and SP
- ▶ Adjust "strength" of MP and/or SP
- ▶ Add additives to MP → Interact specifically with analyte, SP, MP
- ▶ MP velocity → Does not alter retention (K or K')

Stationary Phase Most use silica support particles

- ▶ Not great in high pH → Use alumina (high pH resistance)
- ▶ Low pH → SP can come off of support (hydrolyzed) → Use polymer support

Almost always a "pure" SP

- ▶ Not mixed

Use chemistry reaction to anchor SP to wall/surface

Almost always a "monolayer" SP

- ▶ $d_f \sim 0$

Wide range of polarities

- ▶ Use chemical reactions to adjust the SP

Can use chiral SP

- ▶ Separate enantiomers
- ▶ Reasonable environment conditions



By shaylannxd

Not published yet.

Last updated 20th June, 2022.

Page 3 of 13.

Sponsored by [ApolloPad.com](https://apollopad.com)

Everyone has a novel in them. Finish Yours!

<https://apollopad.com>

Stationary Phase (cont)

Silane Reaction Use to anchor/bond silicones to silica surfaces

- ▶ In packing materials (particles)
- ▶ FS capillaries

Use to deactivate silanols

Silanol

- ▶ Very reactive
- ▶ Highly polar
- ▶ Expose on surface of silica

Deactivate silanol

- ▶ Use chloro silane
- ▶ Ex: C18
- ▶ Result in silanization of surface

Residual silanol

- ▶ SP is usually of a different polarity (non-polar)
- ▶ Results in tailing of analytical peak

After reacting surface with SP

- ▶ Use a short chain alkyl
- ▶ Take care of residual silanol

If silanol peaks present

- ▶ Column is old
- ▶ Molecules of SP are desorbs or removed from surface

Silanol Interactions

“Standard” silica → SP support particles

- ▶ Has silanols on surface → Si-OH
- ▶ ~50% of Si-OH are reacted to Si-O-Si-C18

Stationary Phase (cont)

Residual Si-OH

- ▶ When close to a metal in the silica
- ▶ Are “acidic” and deprotonate easily leaving a Si-O⁻ on the surface

Act as ion-exchange sites for basic analyte

- ▶ Reverse-Phase MP is not suited to ion-exchange separation
- ▶ Very poor peaks are obtained for basic analyte → Tailing peak

Alternative options

- ▶ Use a high purity silica column → Less acidic silanol
- ▶ Purchase a deactivated silica column

Affects only basic compounds

- ▶ Neutral and acidic compounds does not show tailing

Particle and Surface Area

Terms that dominate:

- ▶ Overall plate height
- ▶ Overall plate number
- ▶ d_p → VD eq.: A and $C_m U$

Spherical Particle

- ▶ Surface area/Volume scales with $1/d_p$
- ▶ ↑ A/V = ↓ Retention
- ▶ More SP packed = ↑ Retention (K') = ↑

Resolving Power (R')

Size

- ▶ Nearly all SP are μm scale silica particles
- ▶ Impacts VD equations



By shaylannxd

Not published yet.

Last updated 20th June, 2022.

Page 4 of 13.

Sponsored by **ApolloPad.com**

Everyone has a novel in them. Finish Yours!

<https://apollopad.com>

Stationary Phase (cont)

Porous Particles

- ▶ ↑ Surface area per particle
- ▶ ↑ Amount of SP inside column
- ▶ ↑ Retention and sample capacity = Better R'
- ▶ Smaller the pore → the larger the surface area/g of support

Diffusional Trap of small pores

- ▶ Loss of analyte → Tailing peaks
- ▶ Large MW analyte go into small pores → Never gets eluted

Separation molecules

- ▶ Small molecules ~ 80Å
- ▶ Large proteins ~120-300Å

Normal Phase (NP)

Developed initially

- ▶ Used raw silica as SP → Polar silanol

Separation based

- ▶ Polar-polar interactions with silanol
- ▶ Non-polar elute earlier
- ▶ Polar analyte elute later

In general

- ▶ MP is opposite polarity to SP
- ▶ Works well for polar species only

Reverse Phase (RP)

Use non-polar SP → Silane reaction

- ▶ Use polar MP (Water based)

Separation based

- ▶ Nonpolar-Nonpolar interactions
- ▶ Reverse separation of normal phase
- ▶ Polar species elute first
- ▶ Non-polar species elute later

Stationary Phase (cont)

More popular

- ▶ Organic solvents used for MP → Expensive/dangerous
- ▶ Most analytes are made our of biological origin → Soluble in water-based MP

Controlling Retention

SP Retention depends on

- ▶ **Polarities**
 - ▶ Mass of SP
 - ▶ Type of SP

Mass

- ▶ Control by chain length
- ▶ Density of SP on silica → % of silanol reacted
- ▶ Surface area → Porosity

Selectivity depends on

- ▶ Type of SP
- ▶ Chain length
- ▶ Linker type/length

- ▶ **Chain Length**
 - ▶ ↑ Chain length and/or % organic (carbon) load of SP = ↑ Retention (K')

Example

- ▶ C4 → C8 chain
- ▶ Double chain length → Double volume of SP → Double retention

If within the same type of SP

- ▶ Ex: Alkyl chains
- ▶ No significant changes in selectivity
- ▶ Only shrinking/expanding the chromatogram about t_m
 - ▶ Shrink c-gram = reduce carbon
 - ▶ Stretch c-gram = Increase carbon load



By shaylannxd

cheatography.com/shaylannxd/

Not published yet.

Last updated 20th June, 2022.

Page 5 of 13.

Sponsored by [ApolloPad.com](https://apollopad.com)

Everyone has a novel in them. Finish

Yours!

<https://apollopad.com>

Controlling Retention (cont)

If analyte is not retained
 ▶ Increase % carbon of SP/chain length

If analyte is excessively retained
 ▶ Reduce % carbon

Significant changes in selectivity → Resolving power

- ▶ By changing the type of SP
- ▶ Overall retention should remain roughly constant
- ▶ Keeping the % carbon constant

Effect of MP Strength

MP plays an active role in retention
 ▶ Distribution constant

Common solvents can be sorted according to their polarity

Polarity of MP
 ▶ Main factors of controlling K → K'

When changing MP strength
 ▶ Can calculate the retention under new MP

Equation (Only for RP):
 ▶ $K'_{new}/K'_{old} = 10((P'_{new}-P'_{old})/2)$

Equation (Polarity) $P'_{MP} = \text{Weighted Polarity} = \%A * P_{IA} + \%B * P_{IB}$
 ▶ PI= Polarity

Equation (Only for NP):
 ▶ $K'_{old}/K'_{new} = 10((P'_{old}-P'_{new})/2)$

Example

1. Look at K' of first and last peak
 ▶ $K'_{old} = (\text{last peak} - \text{first peak})/\text{first peak} \rightarrow K' = (2.8-1.8)/1.8 = 0.5$
 ▶ $K'_{new} \rightarrow$ Want it at 10
2. Replace terms in equation
 ▶ $K'_{new}/K'_{old} = 10/0.5 = 20$

Controlling Retention (cont)

3. Old MP polarity
 ▶ If old MP is 20% water/80% Acetonitrile
 ▶ $P'_{old} = (0.2)(10.2) + (0.8)(5.8) = 6.68$

4. New MP polarity
 ▶ $P'_{new} = 2\log(K'_{new}/K'_{old}) + P'_{old}$
 ▶ $P'_{new} = 2\log(20) + 6.68 = 9.28$

5. Solve new MP components
 ▶ $P'_{new} = (x)(10.2) + (1-x)(5.8) = 9.98$
 ▶ 79.1% water and 20.9% ACN

Rule of 3

- ▶ Tool to check/estimate results (not used in calculations)
- ▶ Change of 20% water ~ 3x change in K'

MP

Dynamically adjust MP

Gradient

Some sample contain wide range of analytes

- ▶ Low or high retention
- ▶ no single MP that will elute them all in a satisfactory range of k'

MP gradient

- ▶ MP strength is initially "weak" → Analyte well retained
- ▶ Those with low retention → Elute at reasonable K'
- ▶ Strengthen MP over the course of separation
- ▶ Strongly retained species can be eluted → At a reasonable K' and R'

MP ≠ constant → Changing strength

- ▶ K and K' ≠ constant
- ▶ Can no longer be predicted



By shaylannxd

Not published yet.
 Last updated 20th June, 2022.
 Page 6 of 13.

Sponsored by [ApolloPad.com](https://apollopad.com)
 Everyone has a novel in them. Finish Yours!

Controlling Retention (cont)

MP Selectivity Alter selectivity (α) by changing type of solvent

Resolving power

- ▶ $R = (\alpha - 1)(K'/(1 + K'))(\sqrt{N}/4)$
- ▶ Sensitive to selectivity → Critical pairs in peaks

Selectivity

- ▶ Depends on nature of MP
- ▶ Change selectivity = change type of MP

Try and keep

- ▶ $P_{old} \sim P_{new}$
- ▶ Overall retention is roughly the same
- ▶ Selectivity of peaks change

Selectivity changes cannot be predicted

Chiral Separation

Basic Important to bioanalyses and pharmaceutical separation

Theory

Separation of chiral species

- ▶ Enantio selective

Chiral SP or chiral additives to MP → Separation of enantiomers

Possible to separate structural isomers → Strength of interaction changes as a function of isomer

Ion-Exchange

Ionic Species Small "hard" ions

- ▶ Inorganic ions
- ▶ Cannot use ion pairing
- ▶ Ions can interact with appropriate SP → Ionic SP

Ion-Exchange (cont)

Basic Ion-exchange

Theory

- ▶ Equilibrium-based separation
- ▶ Discreet soption and displacement process
- ▶ Carry throughout column

Column

- ▶ Does not use silica particles
- ▶ Use polymer resin
- ▶ Attach with a strong anion or cation

Example Using a strong anion SP

- ▶ Cation exchange column

1. SP sulfonic acid (IEX resin) → Anion surface particle (SO_3^-)

- ▶ Wash column with acid solution → Cations (H^+)

Cations (H^+)

- ▶ Sulfonic acid is protonated (SO_3H)

2. Inject sample with cation analytes → Metal ion (M^{2+})

- ▶ Metal ions interact with SP
- ▶ Metal ions displace some of H^+ from resin

3. Unbind analyte from SP

▶ Introduce a higher concentration of protons behind analytes → MP gradient

- ▶ H^+ displace weakly bound analytes the move onto strongly bound analyte (cation) → Exchange process
- ▶ Cation analytes is displace off of surface and solubilize in MP



By shaylannxd

Not published yet.

Last updated 20th June, 2022.

Page 7 of 13.

Sponsored by [ApolloPad.com](https://apollopad.com)

Everyone has a novel in them. Finish Yours!

<https://apollopad.com>

Ion-Exchange (cont)

- Analytes move down the column in strength of MP
 - ▶ Each type of analyte elute as a peak
 - ▶ MP ahead of each analyte is too weak → Bound
 - ▶ MP behind each analyte is too strong → Fully displace

Equilibrium Constant

If it behaves like an equilibrium → There is an equilibrium constant

- ▶ Expect to behave like an LC → Produce peaks

Ion-exchange equilibrium constant would behave like a distribution constant

- ▶ Obtain similar result of chromatogram peaks
- ▶ MP controls retention

Equation

$$K_{iex} = \frac{[exchange\&analyte]_s}{[analyte]_m} = C_s/C_m$$

Optimization

Process of Separation

- Carry out initial separation
 - ▶ Choose a strong MP
 - ▶ Ensure everything is eluted and fast separation
- Adjust MP strength
 - ▶ Retention of last peak is within the right region
 - ▶ Depending on the complexity of sample
 - ▶ Simple sample → $K' \sim 10$
 - ▶ Complex sample → $K' \sim 20$
 - ▶ Do calculations for an estimate adjusting needed

Optimization (cont)

- Examine if peaks are within the acceptable region
 - ▶ Examine if all analytes are well resolved
- Consider if a gradient is required
 - ▶ Presence of large area of empty baseline

- If needed
 - ▶ Switch MP type to alter selectivity
 - ▶ Gradient to reach acceptable retention and resolution

- Consider using additives in MP
 - ▶ Help alter selectivity

If MP type/mix strength does not achieve required separation

- ▶ Change SP type
- ▶ May consider type of separation

Summary of MP Effects

Very powerful tool → Versatile

- ▶ Control retention and selectivity

Directly affects distribution constant

↑ MP strength = ↓ K'

MP "strength" is polarity

- ▶ Effect are opposite in RP vs NP
- ▶ RP → Non-polar solvent (organic) = Stronger solvent
- ▶ NP → Polar solvent = Stronger solvent

Ramped MP → Gradient

- ▶ Helps dynamically adjust K'

Useful to make separation less intuitive → $R \downarrow = K' \downarrow$

- ▶ Gradient runs → $R \uparrow = K' \downarrow$



By shaylannxd

cheatography.com/shaylannxd/

Not published yet.

Last updated 20th June, 2022.

Page 8 of 13.

Sponsored by [ApolloPad.com](https://apollopad.com)

Everyone has a novel in them. Finish Yours!

<https://apollopad.com>

LC Detectors

- Ideal Detectors**
- High sensitivity
 - ▶ Steep slope
 - High stability
 - ▶ Minimal drift
 - ▶ Minimal noise on baseline
 - Very low DL
 - Long LDR
 - Can accept MP over wide range
 - ▶ Need reference to null out MP gradients
 - Fast response
 - ▶ Independent of MP
 - Easy to use, maintain and repair
 - Inexpensive
 - Selective/universal
 - ▶ Can be either depending on properties
 - Non destructive
 - ▶ Can collect fractions

- 1λ: UV- Vis Detector**
- Volume
 - ▶ ~1-10 uL (very small)
 - ➔ If V is too large, the signal becomes constant and we see a square shaped peak
 - Pathlengths
 - ▶ ~5-10mm
 - ▶ Longer = better
 - Bigger absorbance for same concentration (beer--lambert law)
 - Window material
 - ▶ Quartz

LC Detectors (cont)

- D2 lamps
 - ▶ Good broad UV source
 - ▶ 185-400nm
 - ▶ Spectrometer to isolate narrow band of wavelength
 - ▶ More simple than FAA
 - ➔ UV does not have to compensate for a flame

- 2 sensors
 - ▶ Sample diode (I)
 - ➔ Intensity coming through the sample
 - ▶ Reference diode (I₀)
 - ➔ Intensity from the light source
 - ▶ Equation
 - ➔ $A = -\log(I/I_0) = -\log(T)$

- Chromatogram
 - ▶ Abs vs time
 - ▶ Use peak area for quantitation

- Many λ: Photodiode array (PDA) Detector**
- Chromatogram
 - ▶ Collect many chromatograms across many wavelength (a spectrum)
 - Sensitivity
 - ▶ Can choose/use chromatogram that provides the greatest sensitivity for each analyte
 - ▶ Find wavelength where analyte has least interference from neighbouring peaks



By shaylannxd

cheatography.com/shaylannxd/

Not published yet.

Last updated 20th June, 2022.

Page 9 of 13.

Sponsored by [ApolloPad.com](https://apollopad.com)

Everyone has a novel in them. Finish

Yours!

<https://apollopad.com>

LC Detectors (cont)

Application

- ▶ Useful to verify which peaks is which when MP is changed

Refractive Index (RI) Detector

- Uses refractive index of analyte compared to MP
- ▶ Snell's law
- The rays will bend if there is a mismatch in refractive indices of the outside and the inside
- ▶ When refractive indices match
- Rays not refracted

Chromatogram

- ▶ If RI match (only MP)
- Full intensity reaches sensor
- ▶ If RI does not match (analyte eluting)
- Reduced intensity reaches sensor
- ▶ Plot signal vs time

Properties

- ▶ Universal
- ▶ Sensitivity
- ~3 orders of magnitude less sensitive than UV
- ▶ Absorbance
- Optically silent
- ▶ Reference flow
- Limited gradient capability

LC Detectors (cont)

Evaporative Light Scatter (ELS) Detector

How it works

- ▶ Uses nebulizer to produce aerosol
- ▶ MP evaporates
- Leaves behind analyte fine crystals
- ▶ Scattering of light (usually laser)
- Only when crystals are present

Analyte

- ▶ Needs to produce crystals
- Very low volatility
- ▶ Can work for non-absorbing analytes
- ▶ Response is nearly uniform for all analytes


Buffers (MP)

- ▶ Must be volatile
- Restricts choices
- Can't use inorganic buffers:
- Leads to buffer salts


Better than RI detector

- ▶ Higher sensitivity
- ▶ Longer LDR

UV-Vis Detector Diagram

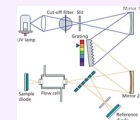

Photodiode Array Chromatogram

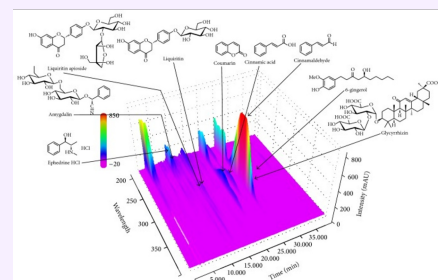

By shaylannxd

cheatography.com/shaylannxd/

Not published yet.

Last updated 20th June, 2022.


Page 10 of 13.

Sponsored by [ApolloPad.com](https://apollopad.com)

Everyone has a novel in them. Finish Yours!

<https://apollopad.com>

Refractive Index Detector Diagram

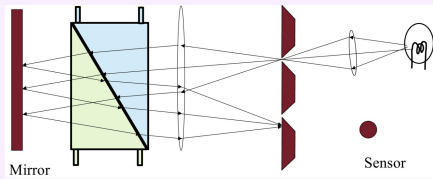

LC-MS (cont)

Mass spectrum

- ▶ Scanning mode → Produce full spectrum → For method development
- ▶ Multiple reaction monitoring (MRM) → Only selected fragments are measured → For quantitation

Better than single quadrupole

- ▶ Lower DL
- Less interferences
- ▶ Longer LDR
- ▶ Allows positive ID of analyte
- ▶ Better selectivity with MRM

Problem

- ▶ Q3 scans across m/z range pretty slowly (1-30 spectra/s)
- ▶ Lowering resolution allows faster scanning → Can't get a full detailed spectrum

QTOF MS

Advantage

- ▶ Can scan 10000 spectra/s
- Many are averaged together to improve quality (better than QQQ)
- ▶ Allows more analytes to be measured simultaneously
- ▶ Higher mass accuracies and resolution
- Permits greater ID power
- ▶ LDR > 5 orders of magnitude

Problem

- ▶ Not as precise as QQQ
- ▶ Expensive

LC-MS

Properties Electropray

- ▶ Sample goes through nebulizer
- ▶ High voltage is applied → Produces charged droplets

Fine metal capillary tube

- ▶ ~0.5-1mm
- ▶ Connected to the outlet
- ▶ Charged with high voltage

Signal

- ▶ MP is pumped
- ▶ Charged droplets are attracted to MS interface
- ▶ Droplets dry down in flight
- charge density ↑ until charge repulsion causes coulombic explosion

Single Mass spectrum

Quadrupole ▶ Simple

MS ▶ MP evaporates away

→ Leaves $[M+H]^+$ ions → no fragments

Problem

- ▶ Difficult for definitive ID
- ▶ Potential m/z overlap

Triple Allows the production of fragments

Quadrupole ▶ Contains Q1, Q2(CID) and Q3

MS

Q1 → Parent ions are selected

Q2(CID) → Collision induced dissociation → Selected ions collide with Ar/He/N₂ (Creates fragments)

Q3 → Fragment ions are filtered/scanned → Detected to produce mass spectrum



By shaylannxd

Not published yet.

Last updated 20th June, 2022.

Page 11 of 13.


Sponsored by **ApolloPad.com**

Everyone has a novel in them. Finish


Yours!

<https://apollopad.com>

Triple Quadrupole Detector Diagram

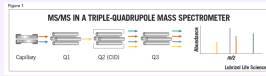

QTOF MS Detector Diagram

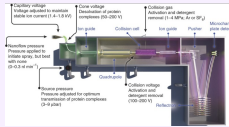

Summary and Applications

Advantage MP plays a critical role in controlling separation

- ▶ Retention → "Strength"
- ▶ Selectivity → "Type"

Wide range of MP available

- ▶ Diverse set of separation conditions
- ▶ Within the same SP and column
- ▶ Allow to quickly try different separation conditions
- ▶ Allow to quickly arrive to a newer optimization separation

No requirements of volatile analyte

- ▶ Needs to be soluble in MP

Wider range of SP available

- ▶ Can choose type
- ▶ Can change particle size
- ▶ Can choose the amount of SP/unit of column

Easy to collect purified analyte

Disadvantage Much lower N compared to GC-FSOT

- ▶ Degrades R and ↑ Overlapping peaks
- ▶ Many LC have low N → 1000-5000

Detector Selective or universal

Comparison

DL

LDR

Summary and Applications (cont)

Cost

- ▶ Purchase
- ▶ Maintenance

Sample capacity

Immune from MP gradients?

Amendable to using IS?

Key Factors if LC is useful

1. Analytes soluble in liquid MP
2. Concentration of analytes are high enough
 - ▶ Can load larger volumes/concentration on columns
 - ▶ Combine with sensitive detectors
3. Does sample require a high R' separation
 - ▶ GC favored over LC
4. Need to recover analyte
 - ▶ LC > GC
5. Slower than GC

Applications

Anti-doping and forensics

- Pharmaceutical
- ▶ Process control
 - ▶ Quality control
 - ▶ R&D
 - ▶ Metabolic
 - ▶ Proteomic

Food and Beverages

- ▶ Vitamins
- ▶ Pesticides
- ▶ Contaminants

Environmental

- ▶ Pesticides
- ▶ Industrial materials



By shaylannxd

Not published yet.
Last updated 20th June, 2022.
Page 12 of 13.

Sponsored by [ApolloPad.com](https://apollopad.com)
Everyone has a novel in them. Finish Yours!
<https://apollopad.com>

Summary and Applications (cont)

R&D

- ▶ Organic synthesis
- ▶ Catalysis

Industrial

- ▶ Feedstock



By **shaylannxd**

cheatography.com/shaylannxd/

Not published yet.

Last updated 20th June, 2022.

Page 13 of 13.

Sponsored by **ApolloPad.com**

Everyone has a novel in them. Finish Yours!

<https://apollopad.com>